

## **Product datasheet for TA150052**

## Mouse monoclonal eGFP antibody, clone OTI5A2

## **Product data:**

**Product Type:** Tag Antibodies

**Product Name:** Mouse monoclonal eGFP antibody, clone OTI5A2

Clone Name: OTI5A2 (formerly 5A2)

Applications: WB

Recommended Dilution: WB:1:500-1:2000

Target: eGFP Host: Mouse

Clonality: Monoclonal

Immunogen: A full length eGFP tagged recombinant protein expressed in HEK293T cell

**Isotype:** lgG<sup>-</sup>

Formulation: PBS (pH 7.3) containing 1%BSA, 50% glycerol and 0.02% sodium azide

**Concentration:** 0.7mg/ml

Conjugation: Unconjugated Storage: Store at -20°C.

**Stability:** Stable for 1 year at -20°C from date of shipment.

**Predicted Protein Size:** 26.9 kDa


**OriGene Technologies, Inc.** 9620 Medical Center Drive, Ste 200

CN: techsupport@origene.cn

Rockville, MD 20850, US Phone: +1-888-267-4436 https://www.origene.com techsupport@origene.com EU: info-de@origene.com



## **Product images:**

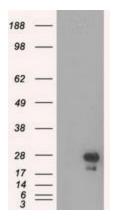

HEK293T cells were transfected with the pCMV6-ENTRY control (left lane) or pCMV6-ENTRY eGFP (right lane) cDNA for 48 hrs and lysed. Equivalent amounts of cell lysates (5 ug per lane) were separated by SDS-PAGE and immuno-blotted with 5A2 anti-eGFP monoclonal antibody (TA150052) (1:2000).

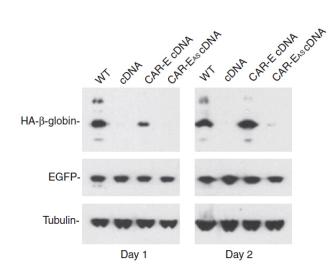


Figure from citation: Western Blot of eGFP-fused protein level by using anti-eGFP antibody in HeLa cells. Dilution: 1:1000 <u>View Citation</u>